

EverFlo Oxygen Concentrator


oxygen simplified



The smaller, lighter EverFlo oxygen concentrator from Respirationics is making a big impression on providers and patients.

Its sleek, ergonomic design makes EverFlo easier to store, move, and carry than conventional concentrators. And it is one of the quietest and most energy efficient devices available today. With no filter for patients to maintain, a durable metal oxygen outlet, and fewer moving parts than other concentrators, EverFlo is designed to deliver consistent, trouble-free operation.



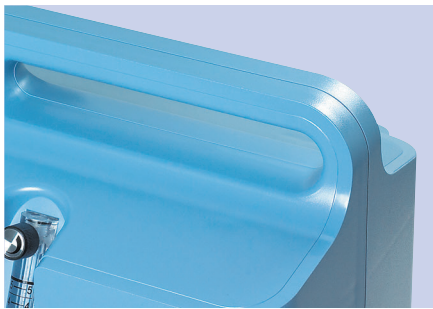

Molded handle that's easy to pick up and transport with one or both hands.



Humidifier bottle holder compatible with all bottle styles.



Internal tubing storage keeps tubing conveniently available for use.



Metal cannula connector for greater durability.



Recessed flow meter reduces the risk of breakage.



Interface panel with on/off switch and alarm indicators.

Product specifications

EverFlo	Standard version	Low power version	IKK version
Power	230V/50Hz	230V/50Hz	230V/50Hz
Liter flow	0.5-5 liters per min	0.5-5 liters per min	0.5-5 liters per min
Sound level	43 dBA	43 dBA	< 40 dBA
Oxygen purity	90% - 96%	87% - 96%	87% - 96%
Weight	14 kg (32 lbs)	14 kg (32 lbs)	15 kg (33 lbs)
Dimensions (HxWxD)	580mm (22.8"), 380mm (15"), 240mm (9.5")	580mm (22.8"), 380mm (15"), 240mm (9.5")	580mm (22.8"), 380mm (15"), 240mm (9.5")
OPI (oxygen percentage indicator) alarm levels	Low oxygen: 82%, very low oxygen: < 70%	Low oxygen: 82%, very low oxygen: < 70%	Low oxygen: 82%, very low oxygen: < 70%
Required maintenance	Compressor intake filter: 2 yrs. Micro disk filter: as needed	Compressor intake filter: 2 yrs. Micro disk filter: as needed	Compressor intake filter: 2 yrs. Micro disk filter: as needed
Power consumption	< 300W	Up to 3 lpm: 285W; < 300W at 5 lpm	< 300W
Part numbers	1020006 Europe 1020009 South Africa 1020011 Italy 1020010 Australia 1020012 Argentina 1020013 South Korea 1020016 Hong Kong 1020017 Switzerland	1020008 UK, Ireland	1020007 Germany

For more information please contact us at:

Air Liquide Healthcare, Stratus House, Unit 1 IDA College Park, Blanchardstown Road North, Blanchardstown, Dublin 15, D15 PEC4.

Tel: +353 (0)1809 1800 Freephone: 1800 24 02 02 (ROI only) 0800 328 5875 (NI only)

Fax: +353 (0)1829 3966 Email: healthie@airliquide.ie Website: www.airliquidehealthcare.ie